


PAJUNK®

Disposable Balloon Systems

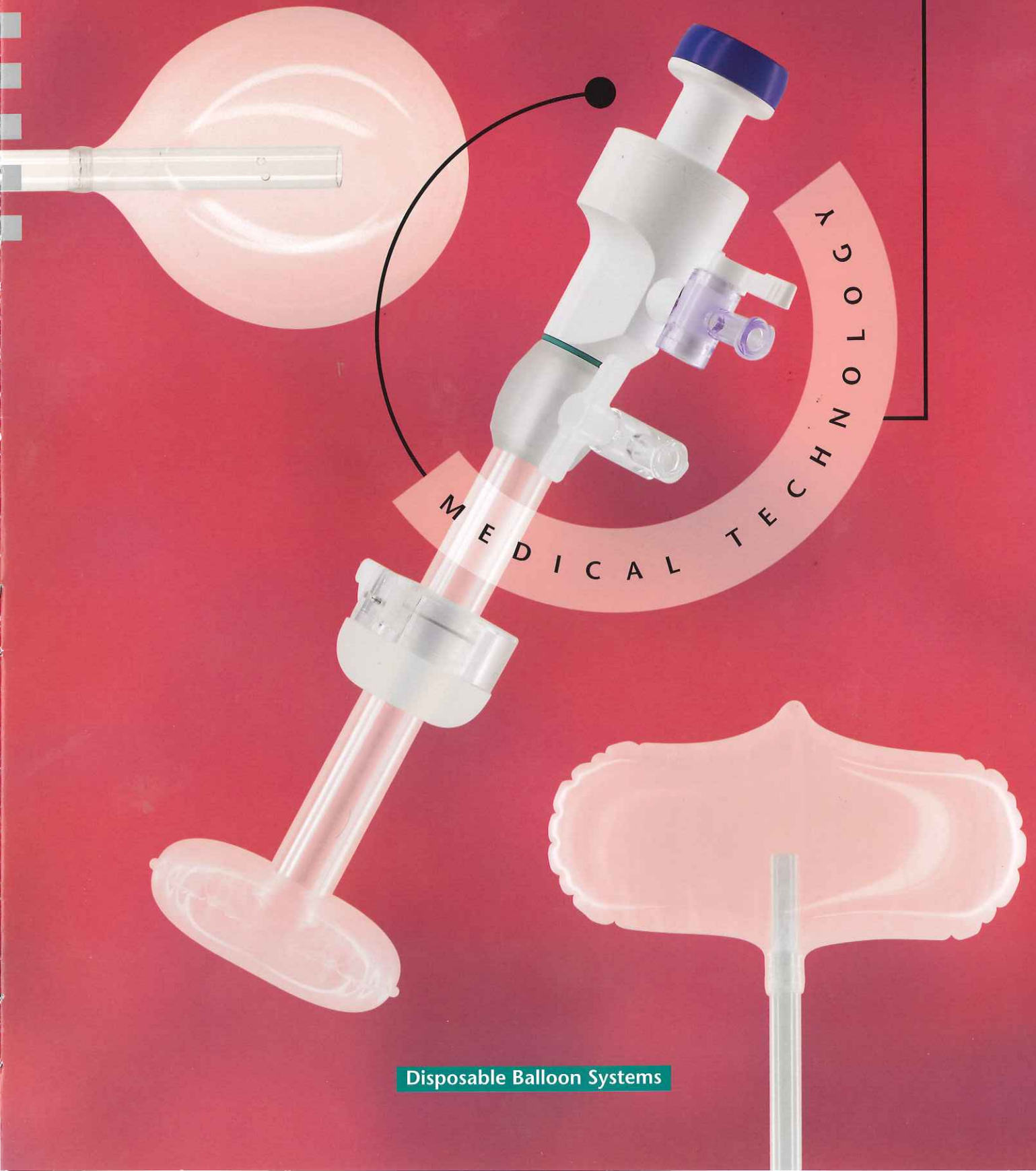

Disposable Balloon Systems

Space Creating Balloons

The space creating PAJUNK® Balloon Systems facilitate optimal, fast and effective dissection of the extraperitoneal space under direct view. Orientation is made easier and three different shaped balloons are available for extra-peritoneal surgeries.




Balloons to Ensure Working Space Access

Structural balloon system

The structural balloon is particularly suitable for extraperitoneal access. The design of the balloon simplifies the separation of the peritoneum from the abdominal wall. Additionally it prevents the extraperitoneal space from collapsing in the case of a loss of gas. Contamination of the lens system with dissected tissue is avoided. The gastight sealing of the incision is formed in the same way as the ring-anchor balloon system.


Ring-anchor balloon system

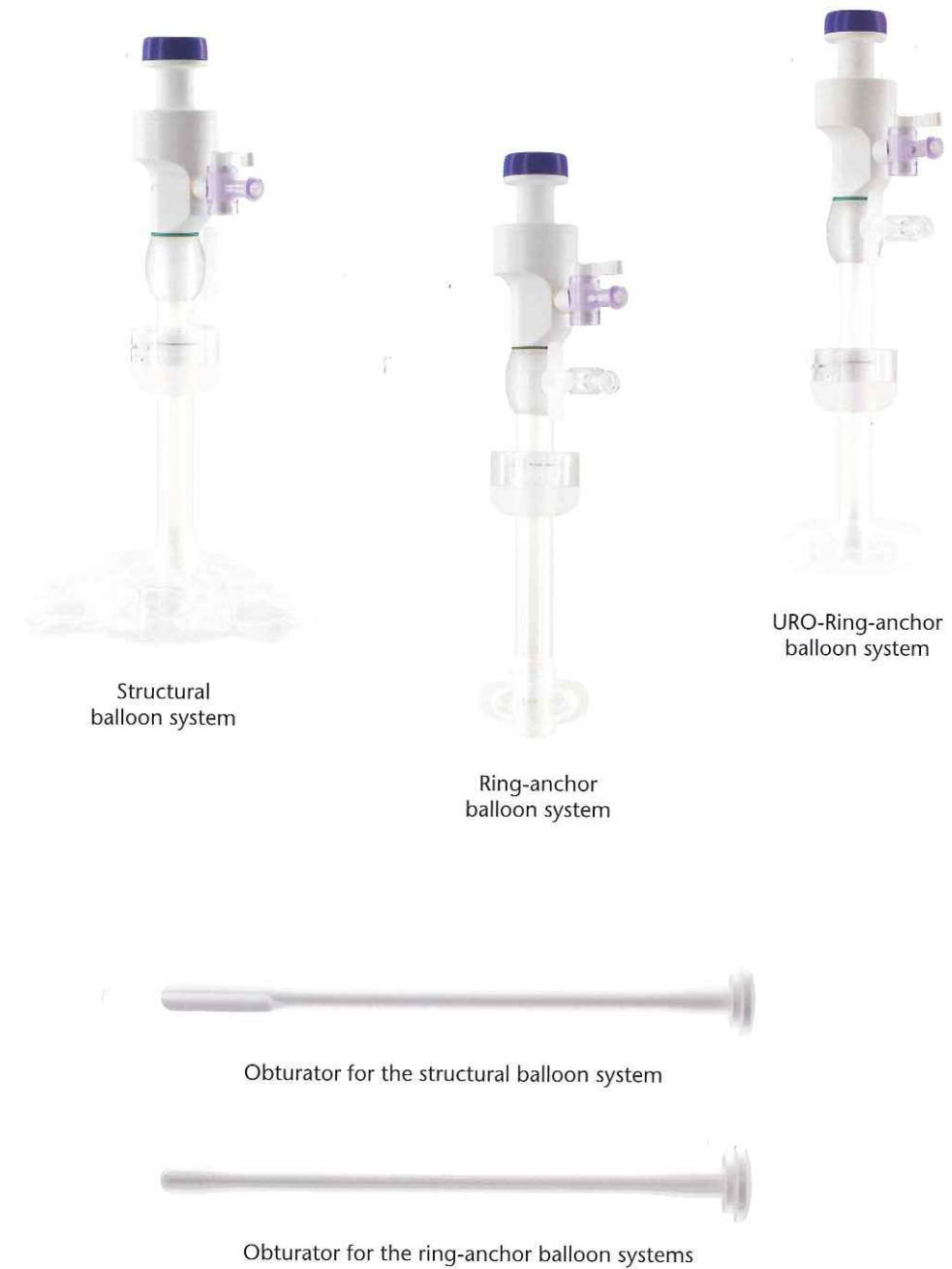


The ring-anchor balloon is also suitable for gastight accessing of the extraperitoneal space. The seal is made by inflating the ring-anchor balloon and then advancing and securing the fastening device. Insufflation is performed via the built-in stop cock.








Balloon Systems-Sets

Space creating balloons

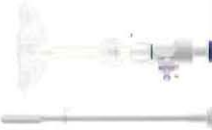

	Product	Size	Order no.	PU*	Ball pump	Syringe 30 ml
	Unilateral dilation balloon Sterile, disposable set with dilation balloon (round)	Ø 11 mm	1284-10-11	5	•	
	Unilateral dilation balloon Extra Large Sterile, disposable set with dilation balloon (round, Extra Large)	Ø 11 mm	1284-60-11	5	•	
	Bilateral dilation balloon Sterile, disposable set with bilateral dilation balloon (oval)	Ø 11 mm	1284-50-11	5	•	

Balloons to ensure working space access

	Product	Size	Order no.	PU*		
	Structural balloon system Sterile, disposable set with structural balloon and instrument channel	Ø 11 mm	1284-20-11	5	•	
	Ring-anchor balloon system Sterile, disposable set with ring-anchor balloon and instrument channel	Ø 11 mm	1284-30-11	5		•
	Ring-anchor balloon system, short Sterile, disposable set with ring-anchor balloon and instrument channel	Ø 11 mm	1284-40-11	5		•
	URO-Ring-anchor balloon system Sterile, disposable set with slim ring-anchor balloon and instrument channel	Ø 11 mm	1284-70-11	5		•
	URO-Ring-anchor balloon system, short Sterile, disposable set with slim ring-anchor balloon and instrument channel	Ø 11 mm	1284-75-11	5		•

PU* = Packaging per Unit